

Siena Care

09000 The Siena nursing bed offers you a very functional and efficient mattress platform. Thanks to its flexible finishing choice, this bed can be perfectly combined with other nursing beds from our gamut.

Carriage

- Legs at a fixed height or optionally adjustable in 3 positions.
- The bed can be optionally equipped with 4 design castors Ø 100 mm, with separate brakes.
- The carriage can be optionally equipped with a central brake mechanism and 4 design castors Ø 125 mm (not height-adjustable).

Mattress platform

- Three- or four- (*) section mattress platform, manually adjustable.
- The backrest is adjustable up to $\pm 70^\circ$, with auto-regression by means of a gas spring of electric actuator.
- The knee rest is adjustable up to $\pm 34^\circ$ by means of a gas spring (in case of 4 sections).
- Available with a mattress platform electrically adjustable (*).
- The legrest or footrest can be manually adjusted.
- The mattress platform consists of rounded steel mattress platform profiles.
- Available with removable solid laminate mattress platform sections of 6 mm (*).
- The backrest and footrest are equipped with 2 synthetic mattress holders.
- Standard equipped with a telescopic bed extension up to + 20 cm.
- The head end is equipped with 2 cases to install an IV rod or lifting pole.
- No visible fastening bolts in the foot panel.
- The bed panels are made of solid oven-dried beech (*).
- Rounded corners to increase the passive safety.

Bed manipulations

- The backrest and knee rest are supported by a gas spring and can be adjusted by means of a lever on both bed sides.
- The electric back- and knee rest (*) is activated by a hand switch.
- The electric finishing can be equipped with a lock box (*), making it possible to switch off the electric functions with 1 single button.

Paneltypes

- Midea (standard)
- Valentino (*)
- Puro (*)
- Felina (*)

Sideguards

- Pollux (*)
- Trix (*)
- Atmosphere (*)

Accessories


We would like to refer to our accessory list.

(*) Option



Dimensions and technical data

- Width: 100,5 cm, Length: 208 cm
- Height of the mattress platform: ± 46 cm, optionally adjustable up to ± 55 cm (in 3 positions).
- Safe working weight: 200 kg
- Weight (without options): 65 kg
- Mattress platform dimensions: 204 x 86 cm
- Mattress dimensions: 195 x 85 cm, thickness: 12 to 15 cm

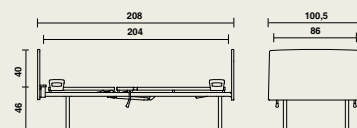

Materials and finishing

- Wood: Massive beech, melaminated chipboard.
- Metal: epoxy-coating RAL 9006 white aluminium.
- Synthetic material.
- Self-lubricating synthetic hinges (= maintenance-free).
- Withstands the common cleaning products.

For more information about the **materials and construction**: please consult our material lists. We would refer to our extensive **colour gamut**.

Options

- 4-section steel mattress platform.
- Mattress platform with removable solid laminate sections.
- Foot extension up to + 20 cm.
- Panel types: Puro, Felina or Valentino.
- Bumpers and/or wall spacer.
- Sideguards.
- Total lock device (ACL).
- Electrical hand switch.
- High/low carriage Vico.
- Removable vacuum-formed synthetic cover cap.
- Wooden surrounds of the mattress platform frame.
- Wooden surrounds for the legs.
- Accessory hooks.
- Blanket rack.
- Legs adjustable in 3 positions.
- 4 design castors Ø 100 mm, with separate brakes.
- 4 design castors Ø 125 mm with central brake (not height-adjustable).

