

Fero Metal recliner chair

06036 The Fero recliner relax chair with lowerable armrests, is a very comfortable chair with specific qualities for the care sector. The reclining mechanism guarantees a long-term seat comfort.

Support

- Firmly welded metal construction.
- All edges and corners have been rounded.
- The legs in front are equipped with synthetic glides.
- Integrated auxiliary wheels at the back legs make it easy to move the chair.

Seat

- Removable seat, made of a melaminated MDF board of 12 mm.
- Upholstered with fire-retardant foam in 2 different density zones to guarantee an excellent seat comfort.
- The rounded seat in front stimulates a fluent blood circulation.
- Comfortable seat height of 45,5 cm.
- The upholstery is finished with imitation leather and decorative stitching.

Backrest

- Made of an ergonomically preformed plywood shell of 12 mm.
- Continuously adjustable up to 24° by means of a gas spring.
- Can be locked in every position, thanks to a lever on both sides.
- Upholstered with fire-retardant foam in 3 different density zones to guarantee an excellent seat comfort.
- The upholstery is finished with imitation leather and decorative stitching.

Lowerable armrests


- The ergonomic front design of the synthetic armrests makes it easy to stand up or sit down.
- Standard equipped with dinner tray attachments.
- Two solid chromium-plated axes to adjust the armrests in height in 4 different positions.
- Locked downwards by means of a spring lever.
- In their lowest position, the armrests are on the same level as the seat.



Dimensions

06036 | Metal recliner chair:

- L 70,5 x D 67-108 x H 115 cm
- Seat height: 45,5 cm
- Seat depth: 51 cm
- Seat width: 53,5 cm
- Height armrests: 66-45,5 cm

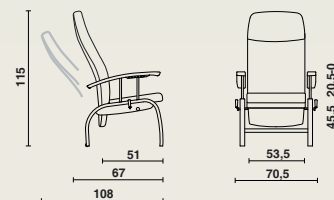

Materials and finishing

- Metal: epoxy coating (light grey), chrome.
- Wood: Plywood, MDF.
- Synthetic.
- Foam: CMHR- and CM-foam (fire-retardant).
- Upholstery: see collection.
- Resistant to the common cleaning products.

For further information about used **materials**, **constructions** and **maintenance**: consult our materialfiles.

Options

- Adjustable and removable head/neck cushion.
- Anatomical reduction cushion.
- Accessory brackets.
- Restraint bracket.
- Lateral headrest.
- Case for IV rod.

Colours

We would like to refer to our extensive colour gamut. Sample card + samples available on request.

Accessories

We would like to refer to our accessorielist.

