Fero Wooden relax chair

07579 The Fero relax chair with castors, is equipped with a relax mechanism that has been especially developed for the care sector.

The <u>Classic mechanism</u> can be easily adjusted by pushing yourself off a little from the chair, and this without using any levers.

Support

- 4 beech laminated legs, connected to each other by means of solid beech support beams.
- 2 beech laminated armrests.
- Glued mortise and tenon joint.
- The sides are firmly connected to each other by means of a solid beech cross beam and a metal frame.
- The ergonomic front design of the armrests makes it easy to stand up or sit down.
- Comfortable rounded armrests of \pm 7,5 cm.
- Standard equipped with a dinner tray attachment.
- All edges and corners have been rounded.
- Four design twin castors Ø 100 mm with separate brakes.

Seat

- Removable seat, made of a melaminated MDF board of 12 mm.
- Upholstered with fire-retardant foam in 2 different density zones to guarantee an optimal seat comfort.
- The rounded seat in front stimulates a fluent blood circulation.
- Comfortable seat height of 47,5 cm.
- The upholstery is finished with imitation leather and decorative stitching.

Backrest

- Made of an ergonomically preformed plywood shell of 12 mm.
- Upholstered with fire-retardant in 3 different density zones to guarantee an excellent seat comfort.
- Equipped with an adjustable and removable head/neck cushion.
- The upholstery is finished with imitation leather and decorative stitching.
- Metal push bar at the back.

Legrest

- Made of a plywood board of 10 mm, upholstered with fire-retardant foam and finished with imitation leather.

Classic mechanism:

- Easy adjustment in 3 positions, i.e. seat position, reading position, and relax position.
- The relax mechanism first raises the legs, after which it is possible to further incline the entire chair.
- Easy to adjust by just pushing yourself off a little, without needing a lever.
- Guaranteed optimal seat comfortin every position: when the chair is reclined, the angle between the backrest
 and the seat remains constant.
- Optimal leg support, thanks to the telescopic legrest.
- When the chair is put in the seat position, the legrest is turned inwards so that you can stand up without any problem.

Colours

We would like to refer to our extensive colour gamut. Sample card + samples available on request.

Accessories


We would like to refer to our accessory list.



פווטופווסוווט

07579 | Wooden relax Fero:

- L 68 x D 82-146 x H 120,5 cm
- Seat height: 47,5 cm
- Seat depth: 51 cm
- Seat width: 53,5 cm
- Height armrests: 64 cm
- Safe working load: 135 kg

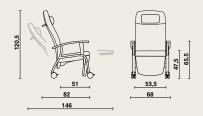

Materials and finishing

- Wood: beech (solid, laminated), Plywood, MDF
- Metal: epoxy coating (light grey).
- Synthetic.
- Fire-retardant foam.
- <u>Upholstery</u>: see collection.
- Resistant to the common cleaning products.

For further information about used materials, constructions and maintenance: consult our materialfiles.

Options

- Seat cushion with visco-elastic foam and a bielastic cover.
- Closed side panels.
- Accessory brackets.
- Restraint bracket.
- Lateral headrest.Case for IV rod.
- Backrest with removable headrest.
- Footrest (height of the seat=52cm).
- Dinner tray support on backrest.

