llana Hospital

10955 The stylish Ilana bedside table is made of very durable materials such as solid laminate and aluminium.

The bedside table is also available in a double-sided version, and can be optionally equipped with a refrigerator and/or dinner tray.

One-sided version (storage cabinet, open compartment, drawer)

- The sides are made of 8 mm solid laminate panels that are assembled between 4 aluminium corner profiles.
- Solid laminate doors.
- Equipped with high-quality self-closing hinges, opening up to 175°.
- Metal design handle.
- An open compartment is accessible from all sides.
- The drawer can be opened on both sides. It is made of steel plate and has 2 solid laminated fronts.
- The drawer glides on 2 synthetic rails.
- The drawer is equipped with a removable synthetic drawer tray.
- Aluminium and steel parts are protected by an epoxy-coating or an anodisation layer.
- Tabletop with rounded solid laminate corners.
- 4 design castors Ø 65 mm.

Optional dinner tray Livia

- Folding
- Easy to incline in both directions.
- The dinner tray has an upright edge.
- The folding part is fixed on a strong aluminium guiding profile that has been finished with an anodisation layer.
- The dinner tray is height-adjustable (between 79 and 112 cm) by means of a gas spring.
- Dimensions of the dinner tray: 63 x 40,5 cm.
- The gas spring is activated by means of a foot pedal.
- The dinner tray is made of fire-retardant ABS.
- The bedside table stands on 5 design castors \emptyset 65 mm.
- The dinner tray can be optionnaly provided with a Click 'n' turn mechanism which allows to move the dinner tray from one side to the other side of the bedside table.

Handles


We would like to refer to our extensive colour gamut. Sample card + samples available on request.

- Open handle

- Closed handle



10955 | 51 x 42,5 x height 86 cm with dinner tray Livia | 67 x 42,5 x height 86 cm

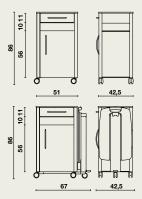

Materials and finishing

- Solid laminate
- Metal: aluminium, steel.
- High-quality self-closing hinges.
- Withstands the common cleaning products.

For more information about the materials, construction and maintenance: please consult our material lists.

Options

- Livia dinner tray.
- Refrigerator.
- Bottle holder.
- Electronic lock on the drawer.
- Cylinder lock on the drawer.
- 2nd drawer in the open compartment.
- Gallery.
- Cases for the telephone holder.
- Removable tray with upright edges for the tabletop.
- 4 bumper castors.
- 4 castors with brakes.
- Click 'n' turn mechanism.
- USB charger and 2 sockets built into the back.



Colours