

DuoSafe sideguards

The DuoSafe sideguards guarantee a flexible fall prevention of the patient. By means of the integrated controls, both the patient and the nursing staff can activate the bed functions.



MOBILITY SUPPORT

The ergonomic side panel offers the patient the necessary support when he wants to get up or walk around the bed. Thanks to the strategic position of the **high-low button**, the mattress platform can be lifted while the patient holds on the sideguards. To limit the risk of getting stuck, there is no gap between the bed and DuoSafe.



COMFORTING DESIGN

Comfort and rest improve the patient's mental condition and speed up the recovery process. Thanks to the openings in the DuoSafe sideguards, the patient can look aside when he is lying in bed. This **visual contact with the environment** is widely considered comforting.



PRIME

Bed panel with a unique handle

The vertical handles make it possible to have a better grip when moving the bed or activating the bed functions. The hands of the nursing staff are protected while moving the bed.



Keypads DuoSafe

KEYPAD - NURSING

All functions within reach

Situated on the outside of the DuoSafe sideguards, the nurse control panel is more tactile than ever before. It is equipped with all functions to put the patient in the desired seating, lying or reading position. Additional functions make it possible for the nursing staff to lock the bed functions and to check the bed position and battery status in the blink of an eye.



KEYPAD - PATIENT STANDARD

Intuitive and within easy reach


The patient panel with integrated background lighting makes it possible to activate the bed functions by means of a touch panel that is within easy reach. The position of the patient buttons makes it possible to activate the bed functions while lying, sitting or getting up.



KEYPAD - PATIENT SMART

Patient empowerment

The intelligent version of the patient panel adds some extra functions to the standard version, such as the activation of room facilities and the possibility to connect to the nurse call system.

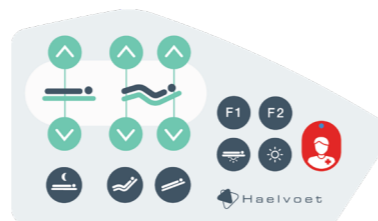

The **background lighting** increases the patient's self-reliance at night. The function buttons illuminate when touched.

Safety first

The keypads of the DuoSafe sideguards are equipped with an indicator that immediately shows you whether the bed is in a safe position or not.

SAFE BED POSITION

Green light



- ✓ The bed is in its lowest position
- ✓ Sideguards in the highest position (option)
- ✓ Central brake is activated (option)




UNSAFE BED POSITION

Red light



- ✗ The bed is **not** in its lowest position
- ✗ Sideguards **not** in the highest position (option)
- ✗ Central brake is **not** activated (option)

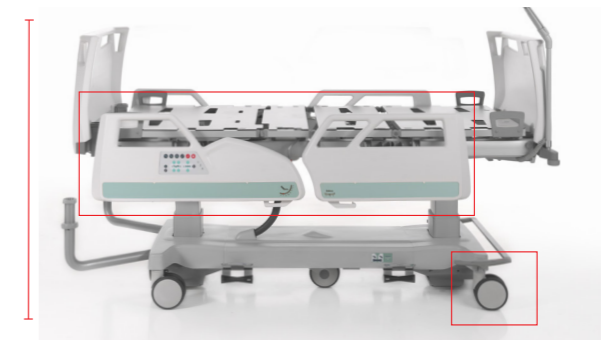

HIDDEN LEVER

with focus on ergonomics

The DuoSafe sideguards can be activated very fluently and easily **with only one hand**. The ergonomic grip has been hidden in the design itself. Thanks to this, **no parts stick out** of the sideguards.

The DuoSafe sideguards are equipped with a special damping mechanism. The movement downwards is slowed down, which increases the patient's comfort and avoids any disturbing noises while resting.




User-friendly design

The DuoSafe sideguards and Prime panels meet the needs of the users. In this respect, a clever design and a durable use of materials is of the utmost importance.

SUPERIOR FUNCTIONALITY

Two **multifunctional hooks** on both sides of the DuoSafe sideguards make it possible to install and use accessories fast and easily. Two integrated **inclination indicators** show the inclination level of both the backrest and the mattress platform. The bed bumpers can be equipped with optional cases for the manual control or accessories.

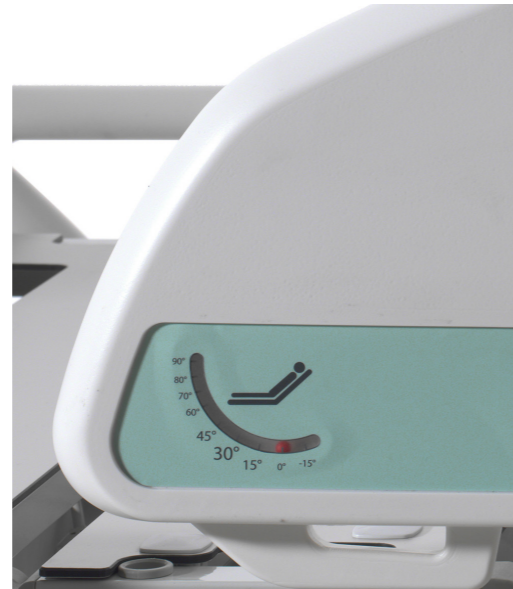

EASY TO CLEAN

Thanks to their clever design and **smooth synthetic surface**, the DuoSafe sideguards can be cleaned very fast and easily. The built-in keypads constitute no obstacle.



CENTRAL LOCKING MECHANISM

The central locking mechanism in the middle of the Prime panel makes it possible to **remove or put back the bed panel fast and easily**.



Bed panels

Our bed structure makes it possible to choose from different bed panel types. Each panel is characterised by its clear lines and keen eye for detail.

		ATMOSPHERE	TRIX	TWIN	DUOSAFE
Type 1: Vela		○	×	×	×
Type 2: Jaro		○	○	○	×
Type 3: Prime		○	×	○	○


○ : Possible × : Not possible



REMOVABLE PANEL

All head and foot panels can be easily removed without needing any tools. The removed **synthetic Prime panel** can stand independently on the ground without any problems.

The decorative strips on the Prime panels and DuoSafe sideguards are standard available in two colours.

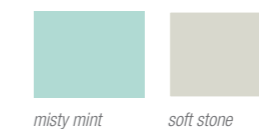

misty mint

soft stone