

Mauro Metal recliner

The Mauro recliner chair with castors, is a very comfortable chair with specific qualities for the care sector.

The reclining mechanism guarantees a long-term seat comfort.

Support

- Firmly welded metal construction.
- 4 Levina design castors Ø 125 mm on precision ball bearings. The 2 castors at the back legs have central brakes with a directional lock.
- All edges and corners have been rounded.
- Extending, synthetic footrest with a step safety and with an antiskid structure (with castors Ø 100 mm and Ø 125 mm).

Seat

- Removable seat, made of a melaminated MDF board of 12 mm.
- Upholstered with fire-retardant foam in 2 different density zones to guarantee an excellent seat comfort.
- The rounded seat in front stimulates a fluent blood circulation.
- Comfortable seat height of 46 cm (up to the footrest).
- The upholstery is finished with imitation leather and decorative stitching.

Backrest

- Made of an ergonomically preformed plywood shell of 12 mm.
- Continuously adjustable up to 24° by means of a gas spring.
- Can be locked in every position, thanks to a lever on both sides.
- Upholstered with fire-retardant foam in 2 different density zones to guarantee an excellent seat comfort.
- The upholstery is finished with imitation leather and decorative stitching.
- Synthetic push bar at the back.

Armrests

- The ergonomic front design of the synthetic armrests makes it easy to stand up or sit down.
- Standard equipped with dinner tray attachments.

Colours

We would like to refer to our extensive colour gamut.
Sample card + samples available on request.

Accessories


We would like to refer to our
accessory list.



Dimensions

20400 | Metal recliner chair:

- L 72 x D 78-155 x H 128 cm
- Seat height up to footrest: 46 cm
- Seat height without footrest: 50 cm
- Seat depth: 51 cm
- Seat width: 53 cm
- Height armrests: 70 cm
- Safe working load: 150 kg

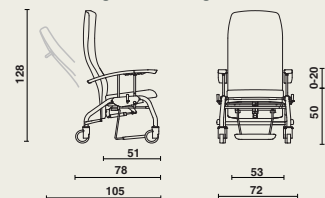

Materials and finishing

- Metal: epoxy coating.
- Wood: Plywood, MDF.
- Synthetic material.
- Fire-retardant foam.
- Upholstery: see collection.
- Resistant to the common cleaning products.

For further information about used **materials, constructions and maintenance**: consult our material files.

Options

- Lock mechanism for the lever.
- Seat cushion with a top layer of 40 mm visco-elastic foam and a bi-elastic cover.
- Accessory hooks.
- Restraint bracket.
- Lateral headrest.
- Case for IV rod.
- Removable headrest.
- Double design brake pedal.
- Dinner tray support on backrest.
- Adjustable and removable head/neck cushion.
- Castors Linea Ø 65 mm with individual brake, without footrest.
- Castors Levina Ø100 mm with individual brake.
- Castors Integral S Ø125 mm, with individual brake on front castors.
- Oxygen bottle holder (Ø 147 mm, H 222 mm)

