

Mauro Metal relax chair

23501 The Mauro relax chair with fixed armrests and castors, is equipped with a relax mechanism that has been especially developed for the care sector.

The Classic mechanism can be easily adjusted by pushing yourself off a little from the chair, and this without using any levers.

Support

- Firmly welded metal construction.
- All edges and corners have been rounded.
- 4 Levina design castors Ø 125 mm on precision ball bearings. The 2 castors at the back legs have central brakes with a directional lock.
- Extending, synthetic footrest with a step safety and with an antiskid structure.

Seat

- Removable seat, made of a melaminated MDF board of 12 mm.
- Upholstered with fire-retardant foam in 2 different density zones to guarantee an excellent seat comfort.
- The rounded seat in front stimulates a fluent blood circulation.
- Comfortable seat height of 46 cm (up to the footrest).
- The upholstery is finished with imitation leather and double decorative stitching.

Backrest

- Made of an ergonomically preformed plywood shell of 12 mm.
- Upholstered with fire-retardant foam in 2 different density zones to guarantee an excellent seat comfort.
- The upholstery is finished with imitation leather and double decorative stitching.
- Synthetic push bar at the back.

Legrest

- Made of a plywood board of 10 mm, upholstered with fire-retardant foam and finished with imitation leather.

Classic mechanism

- Easy adjustment in 3 positions, i.e. seating position, reading position, and relax position.
- The relax mechanism first raises the legs, after which it is possible to further incline the entire chair.
- Easy to adjust by just pushing yourself off a little, without needing a lever.
- Guaranteed optimal seat comfort in every position: when the chair is reclined, the angle between the backrest and the seat remains constant.
- Optimal leg support, thanks to the telescopic legrest.
- When the chair is put in the seat position, the legrest is turned inwards so that you can stand up without any problem.

Armrests

- The ergonomic design of the synthetic armrests makes it easy to stand up or sit down.
- Standard equipped with dinner tray attachments.

Colours

We would like to refer to our extensive colour gamut.
Sample card + samples available on request.

Accessories


We would like to refer to our
accessory list.



Dimensions

23501 | Metal relax chair:

- L 72 x D 78-143 x H 126 cm
- Seat height up to footrest: 46 cm
- Seat height without footrest: 51 cm
- Seat depth: 51 cm
- Seat width: 53 cm
- Height armrests: 72 cm
- Safe working load: 150 kg

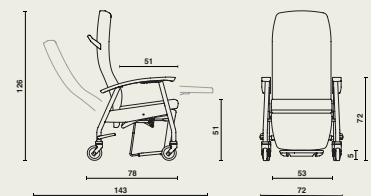

Materials and finishing

- Metal: epoxy coating (light grey)
- Wood: Plywood, MDF.
- Synthetic material.
- Fire-retardant foam.
- Upholstery: see collection.
- Resistant to the common cleaning products.

For further information about used **materials**, **constructions** and **maintenance**: consult our material files.

Options

- Seat cushion with a top layer of 40 mm visco-elastic foam and a bi-elastic cover.
- Accessory hooks.
- Restraint bracket.
- Lateral headrest.
- Case for IV rod Ø 19 or Ø 32 mm.
- Removable headrest.
- Double design brake pedal.
- Anti-static rear castor (with castors Ø125 mm).
- Dinner tray support on backrest.
- Provision for head/neck cushion.
- Adjustable and removable head/neck cushion.
- Castors Ø 65 mm with individual brake, without footrest.
- Castors Levina Ø100 mm with individual brake.
- Castors Integral S Ø125 mm, with individual brake on front castors.
- Raising and sitting assist.
- Oxygen bottle holder (Ø 147 mm, H 222 mm)

