

Forma Hospital

18380 The Forma bedside table offers you an elegant, cost-effective and durable storage unit.

The bedside table is double-sided and equipped with a high drawer and / or dinner tray.

Double-sided version with high drawer (high drawer, open compartment, drawer)

- The sides are made of steel plate.
- Drawers with solid laminate front.
- Metal design handles.
- An open compartment is accessible from both sides.
- The upper drawer is equipped with a removable synthetic drawer tray.
- The high drawer is equipped with a removable tray in synthetic material with place for 3 bottles.
- The metal parts are protected by an epoxy-coating.
- The tabletop is made of solid laminate and has rounded corners.
- 4 design castors Ø 65 mm.


Optional dinner tray Livia

- Folding.
- Easy to incline in both directions.
- The dinner tray has an upright edge.
- The folding part is fixed on a strong aluminium guiding profile that has been finished with an anodisation layer.
- The dinner tray is height-adjustable (between 80 and 110 cm) by means of a gas spring.
- Dimensions of the dinner tray: 63 x 40,5 cm.
- The gas spring is activated by means of a foot pedal.
- The dinner tray is made of ABS.
- The bedside table stands on 5 design castors Ø 65 mm.
- The dinner tray can be optionally provided with a Click 'n' turn mechanism which allows to move the dinner tray from one side to the other side of the bedside table.



Dimension

18380 | with high drawer | 43 x 48 x height 86 cm
with dinner tray Livia | 59 x 48 x height 86 cm

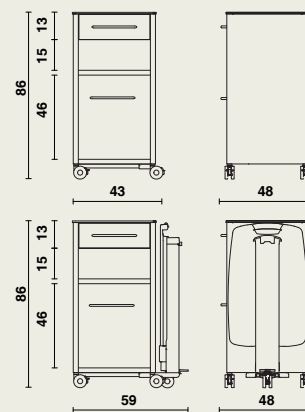

Materials and finishing

- Solid laminate
- Metal: epoxy-coating (light grey or white).
- Withstands the common cleaning products.

For more information about the **materials, construction and maintenance**: please consult our material lists.

Options

- Open compartment finished with solid laminate.
- Livia dinner tray.
- Cylinder lock on top drawer.
- Electronic lock on top drawer.
- Gallery.
- Cases for accessories.
- 4 castors with brakes.
- Click 'n' turn mechanism for Livia dinner tray.

Colours

We would like to refer to our extensive colour gamut.
Sample card + samples available on request.

Handles

- Open handle
- Closed handle

