

# Forma Hospital

**18124** The Forma bedside table offers you an elegant, cost-effective and durable storage unit.

The bedside table is equipped with a high drawer and / or dinner tray.

## Version with high drawer (high drawer, open compartment, drawer)

- The sides and back is made of steel plate.
- Drawers with solid laminate front.
- Metal design handles.
- An open compartment is accessible from on the front side.
- The upper drawer is equipped with a removable synthetic drawer tray.
- The high drawer is equipped with a removable bottle tray (for 3 bottles) in synthetic material.
- The metal parts are protected by an epoxy-coating.
- The tabletop is made of solid laminate and has rounded corners.
- 4 design castors Ø 65 mm.

## Optional dinner tray Livia

- Folding.
- Easy to incline in both directions.
- The dinner tray has an upright edge.
- The folding part is fixed on a strong aluminium guiding profile that has been finished with an anodisation layer.
- The dinner tray is height-adjustable (between 80 and 110 cm) by means of a gas spring.
- Dimensions of the dinner tray: 63 x 40,5 cm.
- The gas spring is activated by means of a foot pedal.
- The dinner tray is made of ABS.
- The bedside table stands on 5 design castors Ø 65 mm.
- The dinner tray can be optionnally provided with a Click 'n' turn mechanism which allows to move the dinner tray from one side to the other side of the bedside table.

## Colours

We would like to refer to our extensive colour gamut.  
Sample card + samples available on request.


## Handles

- Open handle
- Closed handle



## Dimension

**18124 | with high drawer** | 43 x 42,5 x height 86 cm  
**with dinner tray Livia** | 59 x 42,5 x height 86 cm

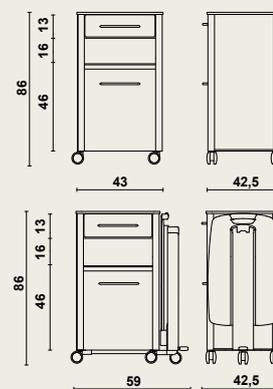

## Materials and finishing

- Solid laminate
- Metal: epoxy-coating (light grey or white).
- Withstands the common cleaning products.

For more information about the **materials, construction and maintenance**: please consult our material lists.

## Options

- Open compartment finished with solid laminate.
- Livia dinner tray.
- Cylinder lock.
- Electronic lock.
- Gallery.
- Cases for accessories.
- 4 castors with brakes.
- Click 'n' turn mechanism.

