


## 2. TWO-PIECE SIDEGUARDS

*Two-piece sideguards not only offer flexibility when needed, but also a lot of advantages when it comes to fall prevention, mobility and maintaining the autonomy of the patients.*

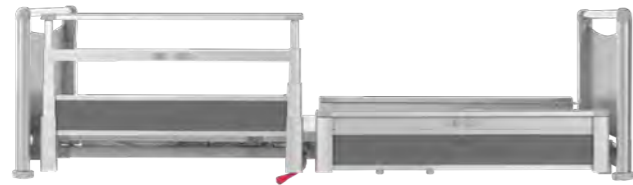

### TWIN

These two-piece sideguards are available as a **single** or **double** version. The double version offers safety over the entire bed length. Twins sideguards can also be used as a raising aid. An integrated safety system prevents the sideguards from being lowered unintentionally.

Max. mattress height Twin sideguards: 18 cm.



### DUOSAFE

DuoSafe panels protect the patient over 3/4 of the bed length. By means of **integrated controls (\*)** both the patient and the nursing staff can activate bed functions and certain optional room facilities.

Max. mattress height DuoSafe sideguards: 19 cm.

**Even when it is extended, the bed complies with safety norm IEC 60601-2:52**, because the gap between the sideguards and the foot section is large enough to prevent somebody from getting stuck (> 318 mm).

(\*): option

### NO RISK OF GETTING STUCK IEC 60601-2-52

All Haelvoet two-piece sideguards comply with safety norm IEC 60601-2:52. This harmonised European norm prescribes the distance between the sideguards and the head section (< 60 mm), between the mattress platform and the sideguards (< 120 mm), between the siderails (< 120 mm), between two-piece sideguards (< 60 mm or > 318 mm) and between the sideguards and the foot section (< 60 mm or > 318 mm).



### TWIN

*Complete safety*

In case **the bed is extended**, a Twin extension piece has to be used to guarantee maximum safety.

