

# Aron+ Junior

**14280** The Aron+ junior bed combines functionality with efficiency for the benefit of both the patient and the nursing staff. Thanks to its lying surface length of 180 cm, this bed is very suitable for intensive or standard care of adolescents.

## Carriage

- Frame with 2 anodised aluminium side beams and 2 steel cross beams.
- High/low adjustment by means of 2 telescopic lifting columns.
- Trendelenburg and Reverse-Trendelenburg adjustment up to  $\pm 16^\circ$ .
- Integral S castors  $\varnothing$  150 mm on precision ball bearings, with central brake mechanism.
- Design pedal to activate the central brake mechanism or directional lock.
- The brake pedal is accessible in all bed positions.
- Finished with a synthetic cover cap.
- Open carriage to improve the accessibility under the bed and the ease of cleaning.

## Mattress platform

- Mattress platform with 4 removable solid laminate mattress platform sections of 6 mm.
- The backrest is electrically adjustable up to  $\pm 70^\circ$  with autoregression.
- The backrest is equipped with an inclination indicator.
- The knee rest is electrically adjustable up to  $\pm 34^\circ$ .
- The footrest is manually adjustable for the Fowler position.
- Electric manipulations of the mattress platform sections are equipped with a clamp safety.
- Manual quick release of the backrest (CPR).
- The perforated, detachable synthetic mattress platform sections can be cleaned very easily.
- Standard equipped with reinforced cutouts for patients restraint straps.
- The backrest and footrest are equipped with integrated mattress holders.
- 4 bed corners with integrated bumper castors.
- 4 bed corners with cases to install an IV rod, lifting pole or accessories.
- Anodised aluminium side beams with integrated rail and accessory hooks.

## Electric bed manipulations

- The electric functions can be activated by means of a hand set or a *mini nurse box ACC* (\*).
- The hand set is wireless and linked to the bed via Bluetooth.
- The bed must be unlocked each time before use, this can be done via the wireless handset or the *mini nurse box ACC* (\*).
- These functions are possible via a key combination on the hand set:
  - Trendelenburg
  - Total lock / unlock (after 3 sec.)
  - CPR
  - personally adjustable transport height
- The hand set is equipped with a pre-programmed sleeping function and button for underbedlight (\*).
- Battery (\*).
- A *mini nurse box* (\*) makes it possible to switch off several functions selectively.
- The software has an automatic lock to prevent accidental use.

### Panel types

- Jaro

### Sideguards

- Atmosphère (\*)

### Accessories


We would like to refer to our accessory list.

(\*) Option



## Dimensions and technical data

- Width: 93 cm, Length: 189 cm
- Height variation of the mattress platform:  $\pm 41$  cm -  $\pm 81$  cm
- Safe working load: 250 kg
- Weight (without options): 125 kg
- Mattress platform dimensions: 175 x 79 cm
- Mattress dimensions: 170 x 78 cm, thickness: 12 to 14 cm

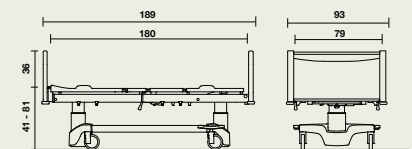

## Materials and finishing

- Metal: epoxy-coating RAL 9006 *white aluminium*, anodised aluminium.
- Synthetic.
- Solid laminate.
- Self-lubricating synthetic hinges (= maintenance-free).
- Withstands the common cleaning products.

For more information about the **materials and construction**: please consult our material lists.

## Options

- Electric protection class I.
- Atmosphere sideguards.
- Mini-nurse box (ACC).
- (Electro-)mechanical quick release (CPR).
- Brake buzzer.
- 5<sup>th</sup> castor.
- Nightlight under the bed (3 LED).
- Battery (lead acid or Li-ion).
- Wall spacer.
- Twin castors  $\varnothing$  150 mm with central brake.
- Inclination indicator.

